

Sleeve print proof — Route B, warm white board

TRUE SIZE · BUILT TO VITALL'S ARTWORK SPEC · 23 AUGUST 2026

Palette variant. Ivory #F4F1EA board, #141414 ink. The research-led default: warm white rather than surgical white, on the grounds that pure white reads clinical in a health frame. **The palette call is still open**, so this proof exists in both; geometry and copy are identical between them.

Route B from **concept-sleeve-fronts-v7-two-routes.html** transposed onto the real dimensions from Vitall's artwork file. Fronts are lifted verbatim from that concept, with one copy change. The back is new: Route B had fronts only, so it has been drawn in the same system and laid out around the compliance window.

THE ONE COPY CHANGE

"Thresholds signed off by Dr Ewa Lindo, GMC-registered GP" is removed from all three fronts. Removed programmatically and asserted at exactly three matches, so it cannot have been missed on one kit. **"ANALYSED IN A UKAS-ACCREDITED UK LABORATORY" stays.**

Worth a decision rather than an assumption: the Ewa attribution still appears on the live site and in approved copy. Taking it off the pack only is fine, but if the intent is that it is not needed anywhere, that is a decision sweep, not a packaging edit.

DIMENSIONS

Flat trim	179.00 × 304.00 mm, from the file's own HiResBoundingBox
Folded face	179.00 × 152.00 mm, lower face rotated 180°
Fold	152.00 mm from trim top — assumed , no crease layer exists in Vitall's file
Bleed	3 mm all round → page is 185 × 310 mm
Die-cut aperture	108.00 × 57.00 mm, 34.09 mm from the left trim edge, vertically centred, corner radius ~3.2 mm (measured off a raster, so a 3.0 or 3.5 mm die spec — confirm)

WHAT THE APERTURE DOES TO ROUTE B

- **It does not touch the fronts.** Route B's marker table, headline and spine are all on the front face, which carries no window (confirmed below).
- **The back is banded around it:** brand and next step above, window clear, regulatory below. There is no room for anything else.
- **The spine flips to the right** on the back in reading orientation, so it wraps the same physical edge as the front's after the fold.
- **No QR on the sleeve.** Per the /activate deprecation the sampling QR belongs on the insert, and the kit is already linked at order creation.

It is a HOLE, not a printed panel (Keith, from a physical kit, 23 August). The board is die-cut; Vitall's kit-details label is stuck to the box underneath and shows through the aperture. That is why their print file leaves the area white and carries no cut layer — the die is a separate tool their kitting house holds. **Matching the existing die exactly is what avoids new tooling cost.**

Confirmed 23 August (Keith): only the back face carries an aperture. The front is clear, so Route B's marker table, headline and spine are unaffected and this proof stands as drawn. The identically sized rectangle on the other face of Vitall's own artwork is their lifestyle photo frame reusing the same dimensions, not a second compliance panel.

WHAT IS STILL OPEN WITH VITALL

ANDROPRIME

ANDROPRIME



AT-HOME BLOOD TEST • TESTOSTERONE

Testosterone

Find out if testosterone is the cause.

TESTOSTERONE	NMOL/L
FREE TESTOSTERONE	PMOL/L
SHBG	NMOL/L
FREE ANDROGEN INDEX	INDEX
ALBUMIN	G/L

ANALYSED IN A UKAS-ACCREDITED UK LABORATORY



CE • UKCA

REF • AP-T01

ANDRO-PRIME.COM


FOLD (assumed) • 152.00

ANDRO-PRIME.COM

CE • UKCA



In-vitro diagnostic for self-collection. Read the enclosed instruction sheet before use.
Distributed by Andro Prime Ltd (Co. 17185839), UK.
Analysed in a UKAS-accredited UK Laboratory.
Store below 25°C. Do not use if the pouch is damaged.

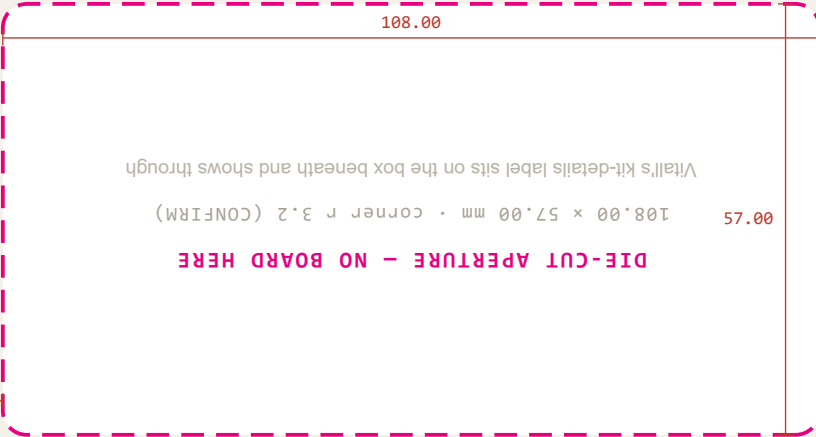

Read the enclosed instruction sheet before taking your sample.

SEE YOUR RESULTS AT ANDRO-PRIME.COM
Your kit is already linked to you. Nothing to register.



ANDROPRIME

ANDROPRIME



AT-HOME BLOOD TEST • ENERGY & RECOVERY

Energy & Recovery

Your levels, not your programme.

VITAMIN D	NMOL/L
ACTIVE B12	PMOL/L
HS-CRP	MG/L
FERRITIN	µG/L

ANALYSED IN A UKAS-ACCREDITED UK LABORATORY



CE • UKCA

REF • AP-E02

ANDRO-PRIME.COM


FOLD (assumed) • 152.00

ANDRO-PRIME.COM

CE • UKCA



In-vitro diagnostic for self-collection. Read the enclosed instruction sheet before use.
Distributed by Andro Prime Ltd (Co. 17185839), UK.
Analysed in a UKAS-accredited UK Laboratory.
Store below 25°C. Do not use if the pouch is damaged.

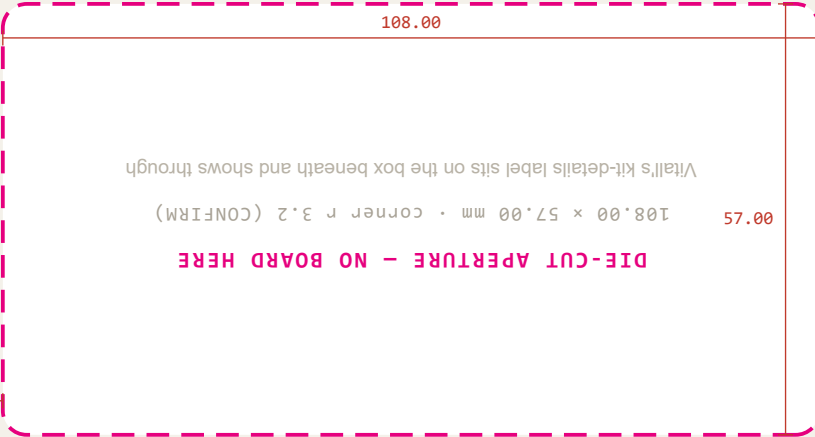

Read the enclosed instruction sheet before taking your sample.

SEE YOUR RESULTS AT ANDRO-PRIME.COM
Your kit is already linked to you. Nothing to register.



ANDROPRIME

ANDROPRIME



AT-HOME BLOOD TEST • NINE MARKERS

Hormone & Recovery

Nine markers. Hormones and recovery in one panel.

TESTOSTERONE	NMOL/L	VITAMIN D	NMOL/L
FREE TESTOSTERONE	PMOL/L	ACTIVE B12	PMOL/L
SHBG	NMOL/L	HS-CRP	MG/L
FREE ANDROGEN INDEX	INDEX	FERRITIN	µG/L
ALBUMIN	G/L		

ANALYSED IN A UKAS-ACCREDITED UK LABORATORY



CE • UKCA

REF • AP-H03

ANDRO-PRIME.COM


FOLD (assumed) • 152.00

ANDRO-PRIME.COM

CE • UKCA



In-vitro diagnostic for self-collection. Read the enclosed instruction sheet before use.
Distributed by Andro Prime Ltd (Co. 17185839), UK.
Analysed in a UKAS-accredited UK Laboratory.
Store below 25°C. Do not use if the pouch is damaged.

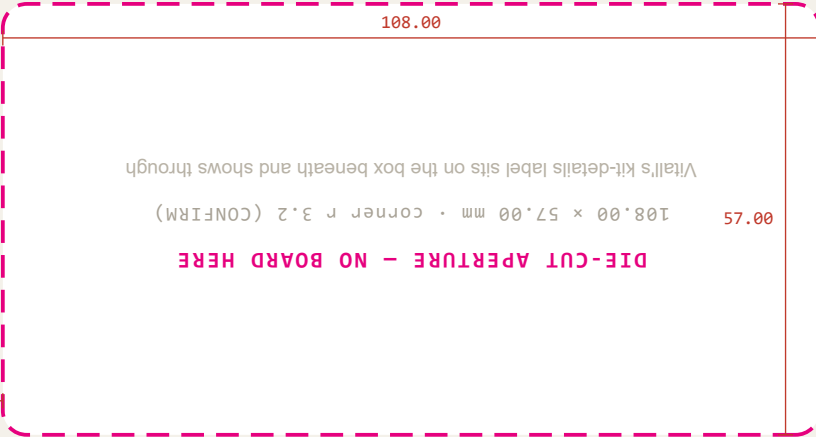

Read the enclosed instruction sheet before taking your sample.

SEE YOUR RESULTS AT ANDRO-PRIME.COM
Your kit is already linked to you. Nothing to register.



Sleeve back — reading orientation

UNIVERSAL ACROSS ALL THREE KITS • TRUE SIZE • MEASURE THIS PAGE TO VERIFY

AP

SEE YOUR RESULTS AT ANDRO-PRIME.COM

Your kit is already linked to you. Nothing to register.

Read the enclosed instruction sheet before taking your sample.

108.00

DIE-CUT APERTURE — NO BOARD HERE

108.00 × 57.00 mm • corner r 3.2 (CONFIRM)

Vital's kit-details label sits on the box beneath and shows through

57.00

34.09

36.91

In-vitro diagnostic for self-collection. Read the enclosed instruction sheet before use.

Distributed by Andro Prime Ltd (Co. 17185839), UK.

Analysed in a UKAS-accredited UK laboratory.

Store below 25°C. Do not use if the pouch is damaged.

IVD

CE • UKCA

ANDRO-PRIME.COM

ANDROPRIME

Band layout, forced by the aperture

Above 0 – 47.50 mm brand + what to do next

Cut 47.50 – 104.50 mm board REMOVED; Vital's label on the box shows through

Below 104.50 – 152.00 mm regulatory fine print

Why the black spine is on the right here

The fold preserves left-to-right, and the back is seen from the other side, so a spine drawn on the right in reading orientation wraps the same physical edge as the front's.